

BODISPA™




THE ALMIGHTY

Cool and Hot massager USER MANUAL



Cold & Hot massager with vibration function. Provides warm or cool treatment for muscle pain relief, stress relief and relaxation. It is great for back, neck, shoulder, arms and legs area.
This massager comes with 2 vibration speeds.

THE PRODUCT DESCRIPTION:

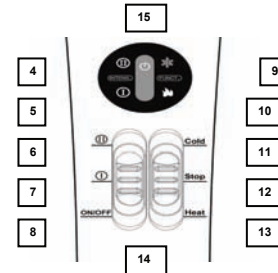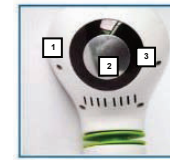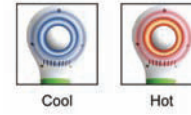

1. Massager Head
2. Hot and Cool treatment disc
3. Signal light
4. Indicating light for power level II
5. Indicating light for power level I
6. power level II
7. Power level I
8. Off function
9. Indicating light for cold treatment
10. Indicating light for hot treatment
11. Cold treatment
12. Stop treatment
13. Hot treatment
14. Switch Button
15. On indicating Light

BODISPA™

BATTERY CHARGE:



Please charge the batteries for approximately 8 hours before initial use.
Make sure the switch in the off position.
Connect the adapter to the AC jack of the unit, and plug in the other end to the electric outlet.
LED indicating light on the adapter will be light up in red while charging. After fully charged, the LED indication light will turn green.

RECHARGING THE BATTERY:

If the power of the massager becomes weak, it means that the battery needs to be recharged.

Note:

Battery life may vary depending on usage conditions

The batteries cannot be overcharged.

HOW TO USE:

Hold the massager and touch on your body where you want to massage.

You can select only the vibration mode or the temperature therapy. You can work both functions at the same time. By pushing the Massage Switch button to select I or II of vibration speed. Or you can also push the Temperature Switch Button to select the hot or cold therapy.

PRECAUTION:



For heat therapy inside THERMOSTAT, will be activated at 40 degree so the blue light signal light of the point no. 3 and the indicating light for cold treatment of point no. 9 will be light up for a while, which is normal.

Do not use when the skin has an injury.

Pregnant women: Do not use this unit!

People who need doctor's special cautions, act upon the doctor's instructions before using this product.

Stop immediately using the unit if unusual symptoms happen when in use

Do not dismantle, repair or modify the unit.

Do not use with wet hands to connect or disconnect the adapter plug from the power outlet.

Do not use a damaged power cord/plug or a defective socket.

Do not use with a helpless person, an infant or a person insensitive to heat

Overlong application at high setting can result in skin burns.

Do not use inside the bath

SPECIFICATION:

Important Instructions, please retain for future use.

Please kindly note that the power for this Cold & Hot massager only to be used with the AC adapter here supplied

Model No:

Input: 100-240V ~ 50/60Hz 10-20VA

Output: 3.5v 800mA

Power consumption: 2.8 watt